



Series 530 Horizontal Sealant Tester

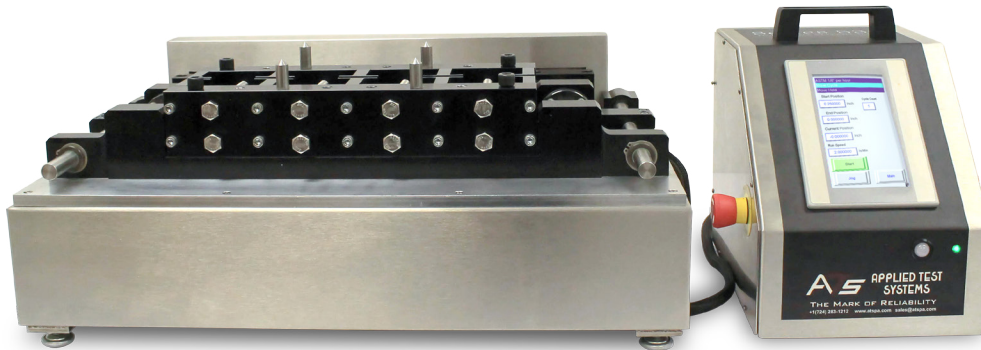

The Series 530 Horizontal Sealant Tester is designed to perform Aymar Type Tests at 1/8 inch per hour at variable speeds of 1/8 – 2 1/2 inches per minute according to ASTM C910, ASTM C736, ASTM C719, ASTM C1305, FEDERAL SS-S200, and FEDERAL TT-S230. We've improved upon our previous design with a touchscreen display, advanced controller options, and user friendly sample fixtures. When operating with an external freezer, users can now easily remove specimens without removing the entire machine. Series 530 fixtures and controllers are interchangeable – one controller can be used on any Series 530 machine, and the fixtures can be easily transferred from one Series 530 to another.

Features

- User friendly sample removal
- Remote connection allows operators to monitor and control their equipment off site
- Touchscreen controller with increased digital range, offering accurate and programmable control
- Increased sample capacity - users can test 8 samples with a standard fixture and up to 16 with the optional additional fixture
- Eight pre-programmed languages available for selection
- Additional sample fixture available for testing to ASTM D 5849

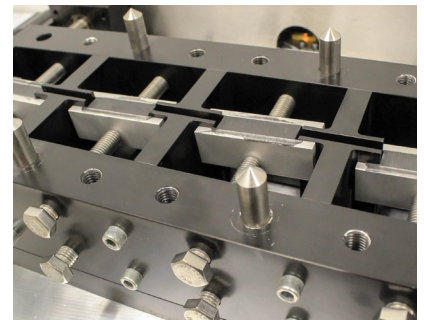

Product Specifications

Number of Cycles	0 - 9,999 cycles
Speed	0 - 5" per minute
Power Requirements	115 VAC, 1 phase, 60 hz or 230 VAC, 1 phase, 50 hz